

Keeping Baby Warm Just Got Easier...




NEW for
Transport



NOW with Integrated Battery
& DC Power for Transport

Conductive Underbody Warming with Flexible
Temperature Settings to Minimize Risk of Cold Stress

Multi-Layer - Multi-Use - Multi-Benefit

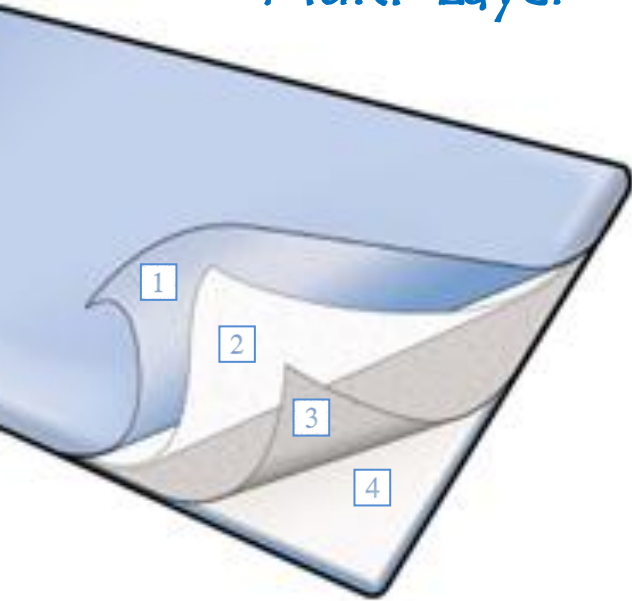

1 Soft, yet Durable Cover

Baby sensitive, Latex-Free. 100% reusable with welded seams for long life. Easy to clean and maintain.

2 Ultra Soothing Layer

Enhanced softness and comfort. Promotes even transfer of warmth without hot / cold spots.

3 Patented Carbon Polymer Thermal Fabric

Delivers continuous, underbody warming at specific, selected temperature setting.

4 Pressure Relief "Memory Foam"

Individualized for maximum support and comfort during this fragile stage of development.

Koala Controller is user-friendly with hi-visibility control buttons for ease of use.

NEW Koala Controller is designed for transport with integrated battery and DC power.

Safety features include thermal cut-off.

Three temperature ranges to choose from:

33 - 37°C / 91 - 99°F

28 - 35°C / 82 - 95°F

32 - 39°C / 90 - 102°F

Settings can be adjusted in 0.5°C / 1.0°F.

Warming Mattress Features

- ▶ Soft, yet durable mattress with patented conductive thermal technology.
- ▶ Delivers active underbody warmth providing a cost efficient alternative to keep baby warm.
- ▶ 100% reusable with no disposable parts. Available in a variety of sizes to fit most bassinets.








Controller Features

- ▶ Compact, energy efficient, lightweight design with hi-visibility controls and adjustable temperature selection available in °C or °F.
- ▶ Provides continuous warming with minimal warm-up time (3-5 minute).
- ▶ Adjustable pole mount for ease of use.
- ▶ Now with integrated battery & DC power .



KOALA^{ΔT} Infant Warmer



PRODUCT DESCRIPTION	MODEL #	DIMENSIONS	WEIGHT
 Standard Warming Mattress	28 - KIW - SM1	24" x 13"	1.3 lbs
 Large Warming Mattress	28 - KIW - LM2	29" x 23"	2.4 lbs
 Warming Blanket	28 - KIW - BL3	18" x 20"	0.9 lbs
 Controller	28 - KIW - CON	6" x 9.5" x 9"	9.0 lbs
 Controller with Integrated Battery and DC Power	28 - KIW - CON - TRANS	11" x 6" x 5"	7.2 lbs

NOVAMEDUSA

4 Westchester Plaza
 Elmsford, NY 10523
 800-425-3535
www.NovaMed-USA.com